

Boehringer Ingelheim	AMTING	PPM Code: 8837903	
		SKU: 145842	
Version: 3	Date: 10.02.2023	Colors	
Product:	<b>Aservo EquiHaler</b>	Black	T
Component Type/Size:	Booklet (343 µg/actuation)	74/14/0/0	<input checked="" type="checkbox"/>
Country(s):	UK	80/0/90/0	<input type="checkbox"/>
Dimensions FLAT (HxWxD):	155 x 180 mm		
Dimensions FOLDED:	n.a.		
Scale/Font Size:	1:1 10 pt.		
Spec/Die/line #:	livret-plié	<b>Finishing / Tech Info</b>	
Manufacturer:	Fareva		
Cancel and Replace:	8837902		
Printer Notes:	n.a.		



# Aservo<sup>®</sup> EquiHaler<sup>®</sup>

343 micrograms/actuation  
inhalation solution for horses  
ciclesonide



Boehringer  
Ingelheim

Boehringer Ingelheim AMT ING		PPM Code: 8837903	
		SKU: 145842	
Version:	3	Date:	10.02.2023
		Colors	
Product:	<b>Aservo EquiHaler</b>	Black	T
Component Type/Size:	Booklet (343 µg/actuation)	74/14/0/0	<input checked="" type="checkbox"/>
Country(s):	UK	80/0/90/0	<input checked="" type="checkbox"/>
Dimensions FLAT (HxWxD):	155 x 180 mm		
Dimensions FOLDED:	n.a.		
Scale/Font Size:	1:1 10 pt.		
Spec/Dieline #:	livret-plié	<b>Finishing /Tech Info</b>	
Manufacturer:	Fareva		
Cancel and Replace:	8837902		
Printer Notes:	n.a.		

Boehringer Ingelheim AMT ING		PPM Code: 8837903	
		SKU: 145842	
Version:	3	Date:	10.02.2023
		Colors	
Product:	<b>Aservo EquiHaler</b>	Black	T
Component Type/Size:	Booklet (343 µg/actuation)	74/14/0/0	<input checked="" type="checkbox"/>
Country(s):	UK	80/0/90/0	<input checked="" type="checkbox"/>
Dimensions FLAT (HxWxD):	155 x 180 mm		
Dimensions FOLDED:	n.a.		
Scale/Font Size:	1:1 10 pt.		
Spec/Dieline #:	livret-plié	<b>Finishing /Tech Info</b>	
Manufacturer:	Fareva		
Cancel and Replace:	8837902		
Printer Notes:	n.a.		



Instruction Video

## Contents

Package leaflet.....	04
Instructions for handling and use of the Aservo® EquiHaler®.....	06
Introduction of the Aservo® EquiHaler®.....	06
Activation of the Aservo® EquiHaler®.....	07
Priming.....	08
Actuation.....	08
Administration.....	10
Fill indicator.....	12
Cleaning of the Aservo® EquiHaler®.....	13
Storage of the Aservo® EquiHaler®.....	14

Boehringer Ingelheim AMT ING		PPM Code: 8837903	
		SKU: 145842	
Version: 3	Date: 10.02.2023	Colors	
Product:	Aservo EquiHaler	Black	T
Component Type/Size:	Booklet (343 µg/actuation)	74/14/0/0	<input checked="" type="checkbox"/>
Country(s):	UK	80/0/90/0	<input checked="" type="checkbox"/>
Dimensions FLAT (HxWxD):	155 x 180 mm		
Dimensions FOLDED:	n.a.		
Scale/Font Size:	1:1 10 pt.		
Spec/Dieline #:	livret-plié	<b>Finishing /Tech Info</b>	
Manufacturer:	Fareva		
Cancel and Replace:	8837902		
Printer Notes:	n.a.		

Boehringer Ingelheim AMT ING		PPM Code: 8837903	
		SKU: 145842	
Version: 3	Date: 10.02.2023	Colors	
Product:	Aservo EquiHaler	Black	T
Component Type/Size:	Booklet (343 µg/actuation)	74/14/0/0	<input checked="" type="checkbox"/>
Country(s):	UK	80/0/90/0	<input checked="" type="checkbox"/>
Dimensions FLAT (HxWxD):	155 x 180 mm		
Dimensions FOLDED:	n.a.		
Scale/Font Size:	1:1 10 pt.		
Spec/Dieline #:	livret-plié	<b>Finishing /Tech Info</b>	
Manufacturer:	Fareva		
Cancel and Replace:	8837902		
Printer Notes:	n.a.		



**343 micrograms/actuation  
inhalation solution for horses  
ciclesonide**

**Marketing authorisation holder**  
Boehringer Ingelheim Vetmedica GmbH  
55216 Ingelheim/Rhein  
Germany

**Manufacturer responsible for batch release**  
Fareva Amboise  
Zone Industrielle  
29 Route des Industries  
37530 Pocé-sur-Cisse  
France

**Active substance and other ingredients**  
Each actuation (ex nostril adapter) contains:

**Active substance:**  
Ciclesonide 343 micrograms

**Excipients:**  
Ethanol 7.9 mg

Clear, colourless to yellowish solution.

**Indication**

For the alleviation of clinical signs of severe equine asthma (formerly known as Recurrent Airway Obstruction - (RAO), Summer Pasture Associated Recurrent Airway Obstruction - (SPA-RAO)).

**Contraindications**

Do not use in known cases of hypersensitivity to the active substance, to corticosteroids or to any of the excipients.

**Adverse reactions**

Mild nasal discharge was commonly observed during safety and clinical studies.

The frequency of adverse reactions is defined using the following convention:

- very common  
(more than 1 in 10 animals treated displaying adverse reactions)

- common  
(more than 1 but less than 10 animals in 100 animals treated)
- uncommon  
(more than 1 but less than 10 animals in 1,000 animals treated)
- rare  
(more than 1 but less than 10 animals in 10,000 animals treated)
- very rare  
(less than 1 animal in 10,000 animals treated, including isolated reports).

If you notice any side effects, even those not already listed in this package leaflet or you think that the medicine has not worked, please inform your veterinary surgeon.

**Target species**

Horse

**Dosage, routes and method of administration**

Inhalation use.

The number of actuations to be administered is the same for all horses. The total treatment duration is 10 days:

**Day 1 to 5:**  
8 actuations (corresponding to 2,744 µg ciclesonide) administered twice daily approximately 12 h apart.

**Day 6 to 10:**  
12 actuations (corresponding to 4,116 µg ciclesonide) administered once daily approximately 24 h apart.

The onset of clinical improvement may take several days. The 10 days treatment schedule should normally be completed. In case of any concerns related to the treatment the responsible veterinarian should be consulted.

The Aservo® EquiHaler® contains sufficient inhalation solution for one horse for the entire treatment duration of the 10 days and an additional amount covering priming and potential losses during administration.

Treatment schedule for use:

Treatment days 1 to 5	Treatment days 6 to 10
8 actuations morning and evening approximately 12 h apart	12 actuations once daily approximately 24 h apart

**Advice on correct administration**

The "Instructions for handling and use of the Aservo® EquiHaler™" is provided in section "Other information" of this leaflet.

**Withdrawal period**

Meat and offal: 18 days.

Not authorised for use in horses producing milk for human consumption.

**Special storage precautions**

Keep out of the sight and reach of children.

This veterinary medicinal product does not require any special storage conditions.

Shelf life after first activation: 12 days.

Do not use this veterinary medicinal product after the expiry date which is stated on the carton and the label after EXP.

**Special warnings**

**Special warnings for each target species**

Special care should be taken when administering the veterinary medicinal product. To ensure an efficacious administration, the breath indicator in the chamber wall of the nostril adapter needs to be observed: when the horse inhales, the membrane of the breath indicator curves inwards. During exhalation, the membrane of the breath indicator curves outwards. The spray should be released at the beginning of inhalation, i.e. when the breath indicator starts curving into the chamber. If the movement of the breath indicator cannot be observed, assure the correct positioning of the nostril adapter. If movement of the breath indicator is still not visible or the movement is too rapid, the product should not be administered. Efficacy of the product has not been established in horses with acute exacerbations (<14 days duration) of clinical signs.

**Special precautions for use in animals**

Safety of the veterinary medicinal product has not been established in horses weighing less than 200 kg body weight, or in foals.

The prescribing veterinarian should assess if the horse has a temperament suitable for a safe and efficacious administration of the Aservo® EquiHaler® in agreement with good veterinary practice. Horses might not adapt to an easy and safe application of the Aservo® EquiHaler® within a couple of days. An alternative treatment should be considered if the horse does not adapt to the treatment with Aservo® EquiHaler®.

The onset of clinical improvement may take several days. The use of concomitant medication (such as bronchodilators) and environmental control may need to be considered in cases of severe clinical signs of respiratory obstruction, at the discretion of the attending veterinarian.

**Special precautions to be taken by the person administering the veterinary medicinal product to animals**

Follow closely the instructions for handling and use of the Aservo® EquiHaler® as provided in the package leaflet section "Other Information".

Acclimatising the horse with a training device, available from the marketing authorisation holder, prior to treatment start has shown to ease the administration of the veterinary medicinal product. Surveillance has shown cases where the horses were non-cooperative and therefore could not be treated according to the product information. In case a horse has a tendency towards defensive behavioural reactions, additional safety precautions could be considered (e.g. employ a second person to handle the horse).

Administration of the product should take place in well ventilated surroundings.

People with known hypersensitivity to ciclesonide or any of the excipients should avoid contact with the veterinary medicinal product.

Inhalative and intranasal corticosteroids may cause rhinitis, nasal discomfort, nosebleed, upper respiratory tract infection and headache. An aerosol filtering mask must be worn during handling and administration. This prevents inadvertent inhalation in case of unintended release of actuations outside the nostril or without the nostril adapter.

The product can cause irritation to the eyes due to its ethanol content. Avoid contact with eyes. In case of accidental eye contact, rinse with large quantities of water.

In case of experiencing an adverse reaction due to accidental inhalation, and in case of eye irritation, seek medical advice and show the package leaflet or the label to the physician.

These precautions should be followed by the person administering the product and persons in close proximity to the horse's head during administration.

The safety of ciclesonide after inhalatory exposure has not been established in pregnant women. In animal studies ciclesonide has been shown to induce malformations in foetuses (cleft palate, skeletal malformations). Pregnant women should therefore not administer the product.

If the Aservo® EquiHaler® is visually damaged it should not be used any more.

It is essential to keep the product out of reach for children.

**Pregnancy and lactation**

The safety of the veterinary medicinal product has not been established during pregnancy or lactation.

Use only according to the benefit/risk assessment by the responsible veterinarian.

The product was shown to be teratogenic following oral administration after high doses in rabbits but not in rats.

**Interaction with other medicinal products and other forms of interaction**

Concomitant use of clenbuterol in a field study in seven horses with severe equine asthma did not indicate any safety concerns.

**Overdose (symptoms, emergency procedures, antidotes)**

After administration of the veterinary medicinal product at up to the 3-fold recommended dose for 3 times the recommended treatment duration no relevant clinical signs could be observed.

**Special precautions for the disposal of unused product or waste materials**

Medicines should not be disposed via wastewater or household waste. Ask your veterinary surgeon how to dispose of medicines no longer required. These measures should help to protect the environment. The cartridge contains residual amount of the product at the end of the course of administration. This should be taken into account at disposal of the used veterinary medicinal product.

**Date on which the package leaflet was last approved**

January 2023

Find more product information by searching for the 'Product Information Database' or 'PID' on www.gov.uk

**Other information**

**Pack size:** One Aservo® EquiHaler® with a nostril adapter and a pre-inserted cartridge. The cartridge contains sufficient inhalation solution for the entire treatment duration (140 treatment actuations) and an additional amount covering priming and potential losses during administration within the 10 day treatment duration. Additionally there is residual solution which cannot be delivered with the required accuracy, and should therefore not be administered. The cartridge cannot be removed from the Aservo® EquiHaler®.

For animal treatment. To be supplied only on veterinary prescription.

POM-V  
Vm 04491/5003

Boehringer Ingelheim AMT ING		PPM Code: 8837903	
		SKU: 145842	
Version: 3	Date: 10.02.2023	Colors	
Product:	Aservo EquiHaler	Black	T
Component Type/Size:	Booklet (343 µg/actuation)	74/14/0/0	<input checked="" type="checkbox"/>
Country(s):	UK	80/0/90/0	<input checked="" type="checkbox"/>
Dimensions FLAT (HxWxD):	155 x 180 mm		
Dimensions FOLDED:	n.a.		
Scale/Font Size:	1:1 10 pt.		
Spec/Dieline #:	livret-plié	Finishing /Tech Info	
Manufacturer:	Fareva		
Cancel and Replace:	8837902		
Printer Notes:	n.a.		

Boehringer Ingelheim AMT ING		PPM Code: 8837903	
		SKU: 145842	
Version: 3	Date: 10.02.2023	Colors	
Product:	Aservo EquiHaler	Black	T
Component Type/Size:	Booklet (343 µg/actuation)	74/14/0/0	<input checked="" type="checkbox"/>
Country(s):	UK	80/0/90/0	<input checked="" type="checkbox"/>
Dimensions FLAT (HxWxD):	155 x 180 mm		
Dimensions FOLDED:	n.a.		
Scale/Font Size:	1:1 10 pt.		
Spec/Dieline #:	livret-plié	Finishing /Tech Info	
Manufacturer:	Fareva		
Cancel and Replace:	8837902		
Printer Notes:	n.a.		

## Instructions for handling and use of the Aservo® EquiHaler®

Please read the following instructions carefully prior to first use of the Aservo® EquiHaler® which can be also found when using the URL [info.equi-haler.com](http://info.equi-haler.com) or the enclosed QR code. The Aservo® EquiHaler® is a product for inhalation for horses. The Aservo® EquiHaler® contains sufficient inhalation solution for one horse for the entire duration of the 10 days treatment and an additional amount covering priming and potential losses during administration.



## Introduction of the Aservo® EquiHaler®

The Aservo® EquiHaler® is for **left hand** use only. While holding the Aservo® EquiHaler® with your left hand, you hold and control your horse with your right hand.



Remove the Aservo® EquiHaler® from the outer carton.



Familiarise yourself with the Aservo® EquiHaler®.

It consists of:

- A** Nostril adapter
- B** Breath indicator
- C** Air inlet
- D** Handle
- E** Prime and release lever
- F** Piercing element with fill indicator

In the following chapters, **activation, priming, administration, cleaning, and storage** of the Aservo® EquiHaler® is described in detail.

## Activation of the Aservo® EquiHaler®



Activation of the Aservo® EquiHaler® has to be performed **only once prior to first use**.



To **activate the Aservo® EquiHaler®**, the piercing element **F** needs to be inserted into the handle of the inhaler without pressing the lever **E**.

Put your right hand under the dark grey piercing element **F** ...



... and push the piercing element **F** up completely into the handle **D** until you will hear a click.

In the final position, the piercing element will disappear entirely in the handle **D** and is no longer visible.

The Aservo® EquiHaler® is now activated, but not ready to use.

An aerosol filtering mask must be worn during handling and administration. This prevents inadvertent inhalation in case of unintended release of actuations outside the nostril or without the nostril adapter.

Boehringer Ingelheim AMT ING		PPM Code: 8837903	
		SKU: 145842	
Version: 3	Date: 10.02.2023	Colors	
Product:	Aservo EquiHaler	Black	T
Component Type/Size:	Booklet (343 µg/actuation)	74/14/0/0	☐
Country(s):	UK	80/0/90/0	☐
Dimensions FLAT (HxWxD):	155 x 180 mm		
Dimensions FOLDED:	n.a.		
Scale/Font Size:	1:1 10 pt.		
Spec/Dieline #:	livret-plié	Finishing /Tech Info	
Manufacturer:	Fareva		
Cancel and Replace:	8837902		
Printer Notes:	n.a.		

Boehringer Ingelheim AMT ING		PPM Code: 8837903	
		SKU: 145842	
Version: 3	Date: 10.02.2023	Colors	
Product:	Aservo EquiHaler	Black	T
Component Type/Size:	Booklet (343 µg/actuation)	74/14/0/0	☐
Country(s):	UK	80/0/90/0	☐
Dimensions FLAT (HxWxD):	155 x 180 mm		
Dimensions FOLDED:	n.a.		
Scale/Font Size:	1:1 10 pt.		
Spec/Dieline #:	livret-plié	Finishing /Tech Info	
Manufacturer:	Fareva		
Cancel and Replace:	8837902		
Printer Notes:	n.a.		

## Priming

Priming is required to ensure accurate initial dosing. Priming is performed **only once** and comprises three (3) actuations (see below). The spray will be fully visible after the third actuation.

When pressing the lever **E** of the Aservo® EquiHaler® for the first time, the lower part of the piercing element with the fill indicator **F** will become visible again.

**Do not push the piercing element back up into the device.**

## Actuation

The Aservo® EquiHaler® is **designed for the left hand use only and for use in the left nostril of the horse only**. While holding and operating the Aservo® EquiHaler® with your left hand, you hold and control your horse with your right hand.

Each **actuation** consists of the following two steps (pictures 5 to 8):



Hold the Aservo® EquiHaler® upright in your left hand.



**Step 1:** Press the prime and release lever **E** until it touches the handle and a click can be heard. Release the lever **E** allowing it to slide back into its starting position.



The display of the fill indicator in the piercing element is partially covered with a red flap.



**Step 2:** Press the prime and release lever **E** again with light pressure only until you hear an audible click. Let the lever slide back into its starting position. The spray will be released subsequently into the nostril adapter **A**. The fill indicator now displays the filling level in %.

Once activated the product should no longer be used after 12 days.

### **Please note:**

Should the piercing element be accidentally pushed completely into the handle again, it will automatically slide into the correct position when the Aservo® EquiHaler® is actuated the next time.



Boehringer Ingelheim AMT ING		PPM Code: 8837903	
		SKU: 145842	
Version: 3	Date: 10.02.2023	Colors	
Product:	Aservo EquiHaler	Black	T
Component Type/Size:	Booklet (343 µg/actuation)	74/14/0/0	<input checked="" type="checkbox"/>
Country(s):	UK	80/0/90/0	<input checked="" type="checkbox"/>
Dimensions FLAT (HxWxD):	155 x 180 mm		
Dimensions FOLDED:	n.a.		
Scale/Font Size:	1:1 10 pt.		
Spec/Dieline #:	livret-plié	Finishing /Tech Info	
Manufacturer:	Fareva		
Cancel and Replace:	8837902		
Printer Notes:	n.a.		

Boehringer Ingelheim AMT ING		PPM Code: 8837903	
		SKU: 145842	
Version: 3	Date: 10.02.2023	Colors	
Product:	Aservo EquiHaler	Black	T
Component Type/Size:	Booklet (343 µg/actuation)	74/14/0/0	<input checked="" type="checkbox"/>
Country(s):	UK	80/0/90/0	<input checked="" type="checkbox"/>
Dimensions FLAT (HxWxD):	155 x 180 mm		
Dimensions FOLDED:	n.a.		
Scale/Font Size:	1:1 10 pt.		
Spec/Dieline #:	livret-plié	Finishing /Tech Info	
Manufacturer:	Fareva		
Cancel and Replace:	8837902		
Printer Notes:	n.a.		

## Administration

The nostril adapter should remain in the nostril during the entire administration of the 8 or 12 actuations. If the nostril adapter slides out of the nostril during administration please re-insert into the nostril again.

The Aservo® EquiHaler® should be administered in a well ventilated area.




Hold the Aservo® EquiHaler® in your **left** hand. Make sure that the air inlet **C** is not obstructed.

Stand on the **left** side of the horse so that the horse's head is next to your right shoulder.

Insert the nostril adapter **A** coming from a horizontal position **carefully** into the horse's left nostril, and gently rotate the Aservo® EquiHaler® ...



... into an upright position. Assure that the nostril adapter is inserted in the nasal cavity.

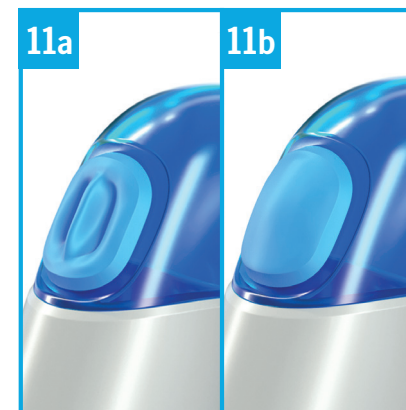

Observe the movement of the breath indicator **B**: When the horse inhales, the membrane of the breath indicator curves inwards (*picture 11a*). When the horse exhales, the membrane of the breath indicator curves outwards (*picture 11b*).

The optimum time for release is at the **beginning** of the horse's inspiration when the breath indicator **B** begins to curve inwards.

**Please note:** In order for the breath indicator to demonstrate when the horse inhales or exhales, the nostril adapter **A** must be correctly placed in the nostril and should fit tightly.

If the movement of the breath indicator cannot be observed, assure the correct positioning of the nostril adapter. If still no movement is visible, the product should not be administered.



Every actuation should be performed following the two steps explained in pictures 6, 7, and 8.

Administer the correct number of actuations as described in section "Dosage, route and method of administration", see above.

Boehringer Ingelheim AMT ING		PPM Code: 8837903	
		SKU: 145842	
Version: 3	Date: 10.02.2023	Colors	
Product:	Aservo EquiHaler	Black	T
Component Type/Size:	Booklet (343 µg/actuation)	74/14/0/0	☐
Country(s):	UK	80/0/90/0	○
Dimensions FLAT (HxWxD):	155 x 180 mm		
Dimensions FOLDED:	n.a.		
Scale/Font Size:	1:1 10 pt.		
Spec/Dieline #:	livret-plié	Finishing /Tech Info	
Manufacturer:	Fareva		
Cancel and Replace:	8837902		
Printer Notes:	n.a.		

Boehringer Ingelheim AMT ING		PPM Code: 8837903	
		SKU: 145842	
Version: 3	Date: 10.02.2023	Colors	
Product:	Aservo EquiHaler	Black	T
Component Type/Size:	Booklet (343 µg/actuation)	74/14/0/0	☐
Country(s):	UK	80/0/90/0	○
Dimensions FLAT (HxWxD):	155 x 180 mm		
Dimensions FOLDED:	n.a.		
Scale/Font Size:	1:1 10 pt.		
Spec/Dieline #:	livret-plié	Finishing /Tech Info	
Manufacturer:	Fareva		
Cancel and Replace:	8837902		
Printer Notes:	n.a.		

## Fill indicator



The **fill indicator** shows the percentage of actuations available in the inhaler. The fill indicator should display 100% prior to first use, i.e. after the Aservo® EquiHaler® is primed.



The display of the fill indicator only moves after several actuations. After administration of the 10 days treatment schedule the display has reached the position **0%**.



The product allows an additional amount of actuations covering potential losses during administration. In this case the display of the fill indicator moves further and stops at the horse head. The inhaler must not be used after the fill indicator has reached the horse head.

## Cleaning of the Aservo® EquiHaler®



After each use and **before cleaning**, check that the fill indicator is blue/white. If it is red, press the prime and release lever **E** until the click is heard. This will assure that you do not accidentally release any spray. To avoid inhalation, hold the inhaler away from your body.



After use, twist and lift nostril adapter **A** from the handle **D**.  
Store the handle in a clean and dry place.



Rinse the nostril adapter **A** **only** in clean running water. Do not use any brushes or cleaning products.


The handle can be carefully wiped with a moist cloth.

The Aservo® EquiHaler® is not suitable for the dishwasher.







 <b>Boehringer Ingelheim</b>	<b>AMTING</b>	<b>PPM Code:</b> 8837903	
		<b>SKU:</b> 145842	
<b>Version:</b> 3	<b>Date:</b> 10.02.2023	<b>Colors</b>	
<b>Product:</b>	<b>Aservo EquiHaler</b>	Black	T
<b>Component Type/Size:</b>	Booklet (343 µg/actuation)	74/14/0/0	<input checked="" type="checkbox"/>
<b>Country(s):</b>	UK	80/0/90/0	<input type="checkbox"/>
<b>Dimensions FLAT (HxWxD):</b>	155 x 180 mm		
<b>Dimensions FOLDED:</b>	n.a.		
<b>Scale/Font Size:</b>	1:1	10 pt.	
<b>Spec/Diefine #:</b>	livret-plié	<b>Finishing / Tech Info</b>	
<b>Manufacturer:</b>	Fareva		
<b>Cancel and Replace:</b>	8837902		
<b>Printer Notes:</b>	n.a.		

Boehringer Ingelheim  
 Vetmedica GmbH  
 55216 Ingelheim/Rhein  
 Germany



8837903